

A polite word from your TA

When writing your reports, please don't lump all the discussion questions into one gigantic incomprehensible paragraph.

Space out the answers instead.

You can even number them if you want.

We sure would love that.

Dr. Anderson makes our lives a living hell, so please be nice to us. You're all we've got left.


Thanks.

CH204

POTIONS CLASS

FALL 2007

PROFESSOR SEVERUS SNAPE

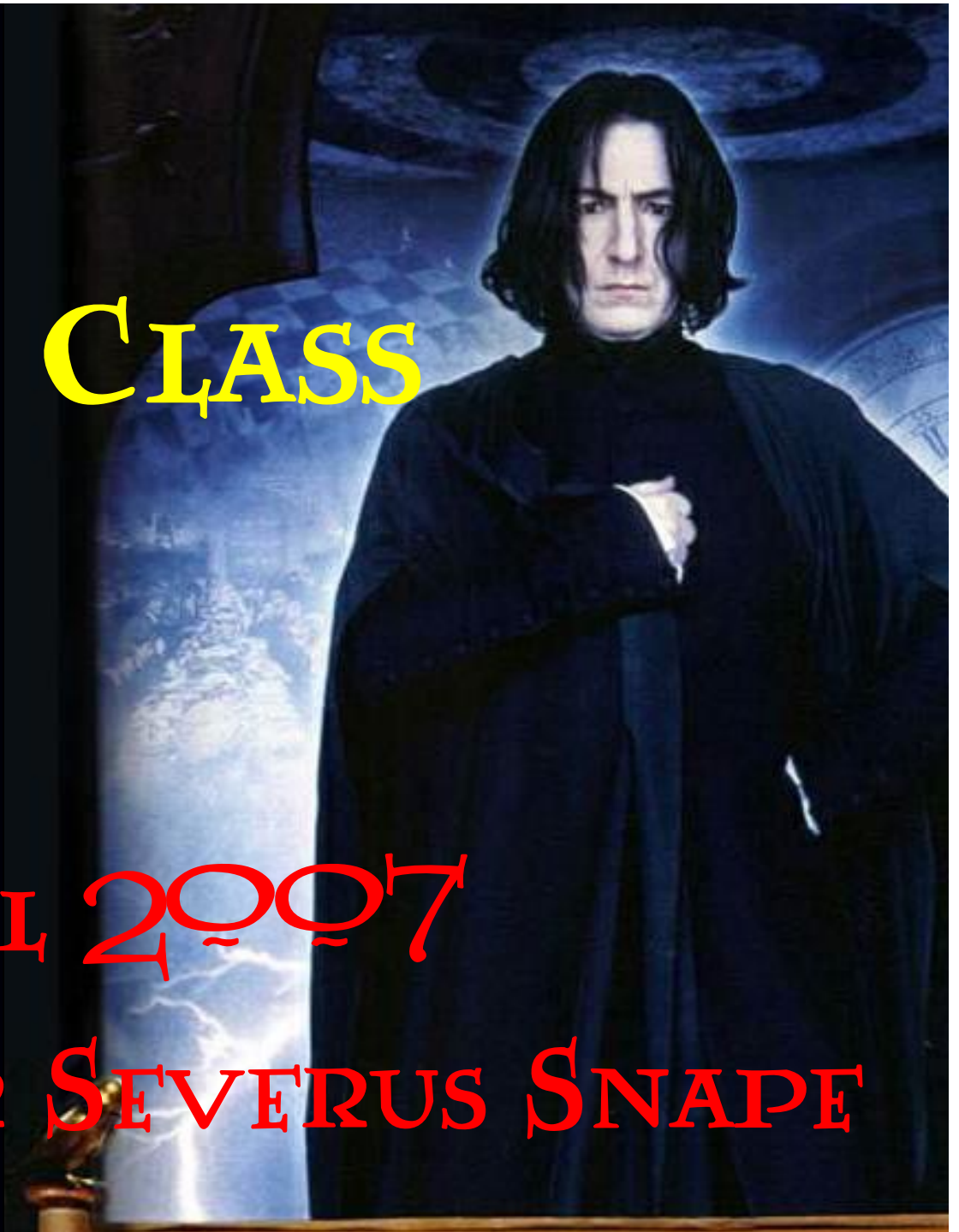

LAST WEEK IN THE POTIONS LABORATORY

SEPARATED MIXTURES BASED ON DIFFERING
PHYSICAL AND CHEMICAL PROPERTIES

USED EXCEL TO CALCULATE AVERAGE, STANDARD
DEVIATION, AND WEIGHT PERCENTS

Q-TEST IN ACTION!

BÜCHNER FILTERING

EXPERIMENT 3
QUALITATIVE CHEMICAL ANALYSIS

“QUANT” VS “QUAL”

QUANTITATIVE – HOW MUCH IS THERE?

QUALITATIVE – WHAT IS IT?

YOU WILL IDENTIFY THE CHEMICAL
IDENTITIES OF 5 UNKNOWN SOLUTIONS
BASED ON HOW THEY REACT
(OR DON'T REACT!)
WITH ONE ANOTHER.

PREVIOUS YEARS

GORGON'S BLOOD

LIQUID GOLD

VITREOUS HUMOR OF A BLIND MULE

2007

ACIDS

BASES

INORGANIC SALTS

TWO-PART LAB

- PART 1: MIX TEN KNOWN SOLUTIONS AND RECORD THE RESULTS OF THE REACTIONS
- PART 2: MIX YOUR FIVE UNKNOWNNS AND COMPARE THE RESULTS WITH WHAT YOU SAW IN PART ONE.

BE EXACT!

**THE MORE ACCURATELY YOU RECORD
YOUR OBSERVATIONS, THE EASIER IT
WILL BE TO IDENTIFY YOUR
UNKNOWN.**

WHAT ARE WE LOOKING FOR?

*EXPLOSIONS

*SUPERNATURAL CREATURES

*RASHES, MUTATIONS,
TRANSFORMATIONS



WHAT ARE WE ^{REALLY} _^ LOOKING FOR?

PRECIPITATES.

(SEE THE SOLUBILITY TABLE IN APPENDIX 2.)

DON'T EXPECT TO SEE ANY ACID-BASE ACTION.

WRITE CHEMICAL EQUATIONS FOR ALL OF THE REACTIONS THAT FORM A PRECIPITATE.

THE KNOWN SOLUTIONS

ACIDS: HCl H₂SO₄ HNO₃

BASES: NaOH Na₂S

SALTS: Ba(NO₃)₂ AgNO₃ K₂CrO₄

Fe(NO₃)₃ Ni(NO₃)₂

ALL SOLUTIONS ARE 0.10 OR 0.20 M.

RIDDLE ME THIS

WHAT DO YOU GET WHEN YOU CROSS
HYDROCHLORIC ACID WITH SILVER NITRATE?

BALANCED CHEMICAL EQUATION



ADD THE PHYSICAL STATES OF EACH COMPOUND




THIS IS CALLED A MOLECULAR EQUATION.

LET'S GET REAL



WRITE AQUEOUS COMPOUNDS AS INDIVIDUAL IONS:

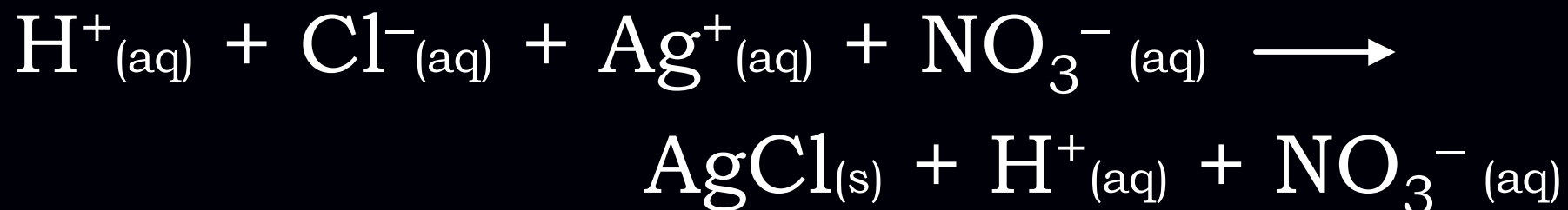


THIS IS A TOTAL IONIC EQUATION.


LOTS OF SPECTATOR IONS.

TIME TO CLEAN HOUSE




CROSS OUT SPECTATOR IONS

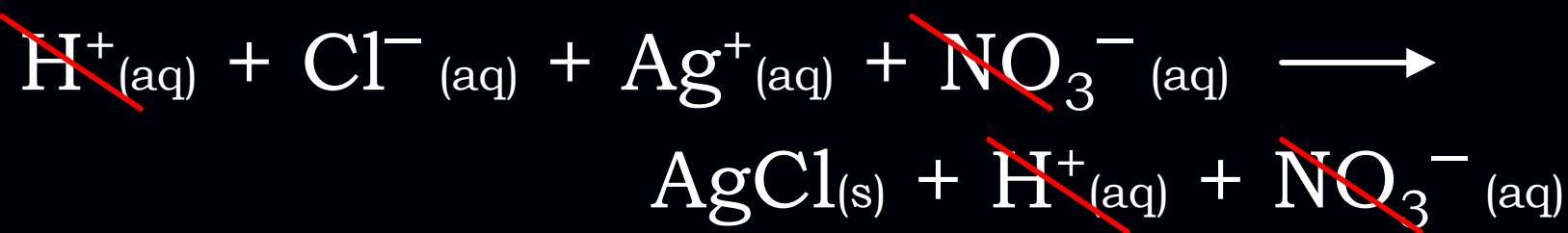


THIS LEAVES US WITH A NET IONIC EQUATION

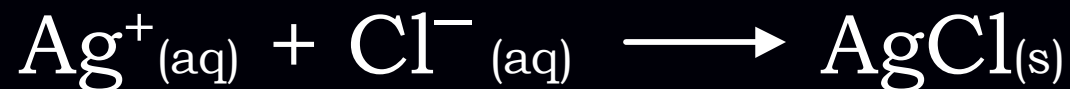

THE NET IONIC EQUATION

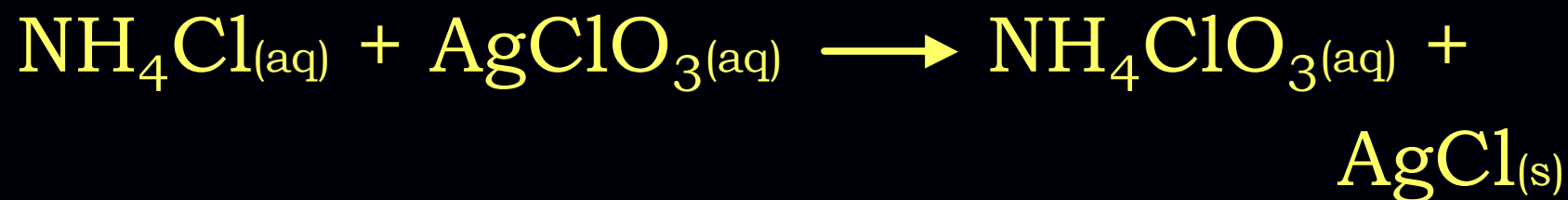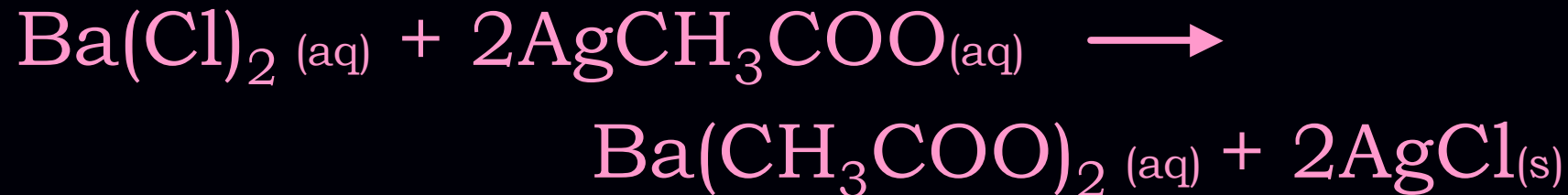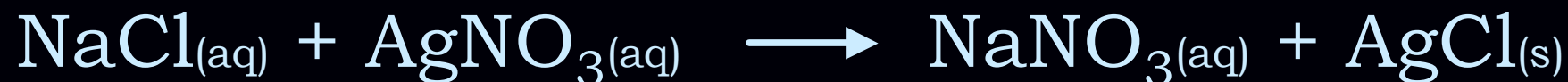

ALL OF THESE REACTIONS HAVE THE SAME NET IONIC EQUATION:

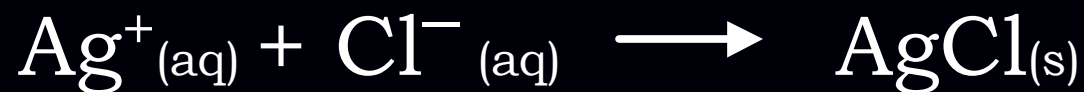

SIMPLE IS GOOD

- THE NET IONIC EQUATION DESCRIBES THE CHEMICAL REACTION THAT OCCURS, AND DOES NOT INCLUDE ANY IONS THAT DO NOT TAKE PART IN THE REACTION, *EVEN THOUGH THOSE IONS ARE PRESENT IN SOLUTION.*
- HOW DO WE KNOW WHICH IONS WILL REACT AND WHICH ONES WON'T?

SOME QUICK SOLUBILITY RULES

- ★ ALL COMPOUNDS CONTAINING ALKALI METALS AND AMMONIUM ION ARE SOLUBLE.




- ★ ALL COMPOUNDS CONTAINING NITRATE, CHLORATE, PERCHLORATE, AND ACETATE ARE SOLUBLE.



SOME QUICK INSOLUBILITY RULES

* ALL COMPOUNDS CONTAINING

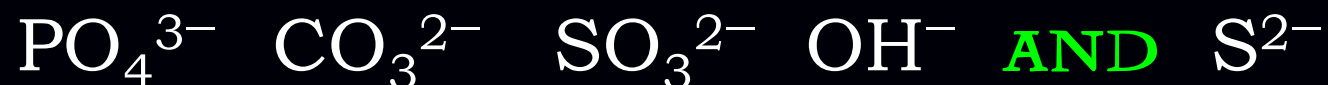

ARE INSOLUBLE EXCEPT THOSE THAT CONTAIN ALKALI METALS, OR NH_4^+ . (EXCEPTION: SOME GROUP II HYDROXIDES ARE SLIGHTLY SOLUBLE.)

* WHEN IN DOUBT, Ag^+ Pb^{2+} AND Hg

COMPOUNDS TEND TO BE INSOLUBLE.

IN THE POTIONS LABORATORY

★ CREATE AN ARRAY OF REACTIONS IN THE MICROWELL PLATE SIMILAR TO THE ONE IN THE LAB MANUAL.

★ USE ONLY 2 DROPS OF EACH REACTANT.

★ *DO NOT* TOUCH THE TIPS OF THE DROPPER BOTTLES TO THE SOLUTIONS IN THE MICROWELL PLATE OR YOU WILL DIE A MOST PAINFUL DEATH.

VILE, HIDEOUS FLUIDS!

EMPTY YOUR USED MICROWELL PLATES
INTO THE DISGUSTING PLASTIC TRAY IN THE HOOD.

RINSE THE PLATES WITH YOUR SQUIRT BOTTLE
AND STACK THE PLATES IN THE HOOD.

LAB REPORT

MOLECULAR EQUATIONS FOR
11 PRECIPITATION REACTIONS.

NET IONIC EQUATIONS FOR
11 PRECIPITATION REACTIONS.

11 + 11 is 22 EQUATIONS ALTOGETHER.

WARNING! DARK MAGIC!

FIVE REACTIONS WILL TURN CLOUDY
EVEN THOUGH NO SOLID SHOULD BE FORMED.



THESE PRECIPITATES ARE DUE TO UNAVOIDABLE TRACE
CONTAMINANTS IN THE SOLUTIONS (POLYSULFIDES IN
 Na_2S AND CARBONATE ION IN NaOH).

SIMPLE LAB, MONSTER WRITE-UP

THE REPORT AND POST-LAB FOR THIS EXPERIMENT
TAKE A LOT OF THOUGHT AND TIME.

ANSWER ALL THE QUESTIONS USING ONLY THE
REAGENTS USED IN THIS EXPERIMENT!

GO TO THE [FREEBIES](#) PAGE OF THE CLASS WEB SITE
FOR A HINT SHEET TO HELP YOU WITH THE
REPORT AND POST-LAB

NEXT WEEK

EXPERIMENT 4: ACID-BASE TITRATION

- MORE INVOLVED THAN THIS LAB.
- REQUIRES MORE PREPARATION THAN OUR PREVIOUS LABS FOR BOTH THE PRELIMINARY WRITE-UP AND THE PRE-LAB.
- **START PREPARING EARLY**

PRE-LAB QUESTION 1: THE ANSWER IS NOT 71!

FINAL EXAM, PART 2

- * YOU WILL NEED A CALCULATOR EVERY WEEK (EXCEPT NEXT WEEK).
- * MAKE SURE YOU KNOW YOUR SECTION NUMBER AND YOUR TA'S NAME!

